



SELF-CATH EASILY WITHOUT RETRACTION

For many users, catheter insertion can be difficult using ordinary closed system catheters, often experiencing the catheter tube slipping back into the bag. With our patented EZ-Advancer® valve, holding the catheter tube in place, users no longer have to worry about retraction during insertion.

Additionally, the MTG EZ-Advancer® closed system catheter makes cathing more hygienic and comfortable, with our soft, flexible introducer tip and fire-polished eyelets for a smooth catheterization experience.

Available in Male, Female, and Pediatric Sizes.

MTG EZ-Advancer®

No-Touch Closed System Catheter

FEATURES & BENEFITS

EZ-Advancer® Valve

Our patented EZ-Advancer® valve secures the catheter tube after each forward advancement, preventing retraction into the bag during catheter insertion.

Designed to Reduce Risk of Recurrent UTIs

The no-touch closed system catheter and pre-lubricated introducer tip are designed to reduce the risk of recurrent urinary tract infections.

Collection Bag Protects Sterility


The self-contained 700ml or 1500ml collection bag with accessible finger holes protects the sterility of the catheter tube and collects urine while measuring output. A sterile sample can be easily collected from the catheter port.

Self-Cath Anywhere with Mini-Pak Size Kits & Pocket-Size Packs


Our Mini-Pak size kits and pocket-size packs are the perfect traveling option and fit discreetly in your pocket or purse. Mini-Pak size kits include a privacy bag for easy and discreet disposal.

Made with Medical Grade Quality Materials

All MTG catheters use medical grade quality materials and are not made with natural rubber latex, BPA or DEHP.

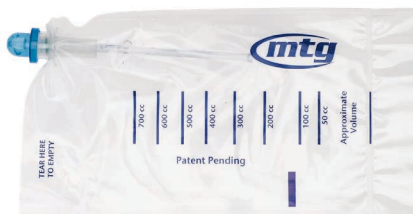

700ml Collection Bag

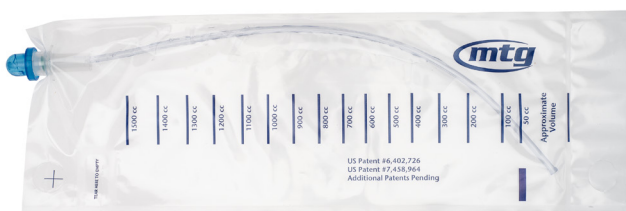

1500ml Collection Bag



MTG catheter kits include everything you need to self-cath. Full kits contain one sterile MTG EZ-Advancer® closed system catheter, pair of latex-free gloves, underpad, gauze pad, Povidone-iodine (PVP) or Benzalkonium chloride (BZK) clear, odor-free skin cleansing prep wipe and privacy bag.

MTG EZ-Advancer® Male No-Touch Closed System Catheters

Item #	FR Size	Length	Description	Antiseptic Wipe
Single Closed System Catheters				
30112	12	16"	MTG EZ-Advancer® No-Touch Closed System Catheter, 1500ml Collection Bag, 12 Fr, 16"	NA
30114	14	16"	MTG EZ-Advancer® No-Touch Closed System Catheter, 1500ml Collection Bag, 14 Fr, 16"	NA
40114	14	16"	MTG EZ-Advancer® No-Touch Closed System Catheter, 1500ml Collection Bag, Pocket-Size, 14 Fr, 16"	NA
Closed System Catheters Kits - Standard Size				
32112*	12	16"	MTG EZ-Advancer® No-Touch Closed System Catheter Kit, 1500ml Collection Bag, 12 Fr, 16"	BZK Swabs
31114*	14	16"	MTG EZ-Advancer® No-Touch Closed System Catheter Kit, 1500ml Collection Bag, 14 Fr, 16"	PVP Swabs
32114*	14	16"	MTG EZ-Advancer® No-Touch Closed System Catheter Kit, 1500ml Collection Bag, 12 Fr, 16"	BZK Swabs
32116*	16	16"	MTG EZ-Advancer® No-Touch Closed System Catheter Kit, 1500ml Collection Bag, 16 Fr, 16"	BZK Swabs
Soft Closed System Catheters Kits - Standard Size				
32212*	12	16"	MTG EZ-Advancer® No-Touch Closed System Catheter Kit, 1500ml Collection Bag, Soft, 12 Fr, 16"	BZK Swabs
32214*	14	16"	MTG EZ-Advancer® No-Touch Closed System Catheter Kit, 1500ml Collection Bag, Soft, 14 Fr, 16"	BZK Swabs
Closed System Catheters Kits - Mini-Pak Size				
42112	12	16"	MTG EZ-Advancer® No-Touch Closed System Catheter Kit, 1500ml Collection Bag, Mini-Pak Size, 12 Fr, 16"	BZK Wipe
42114	14	16"	MTG EZ-Advancer® No-Touch Closed System Catheter Kit, 1500ml Collection Bag, Mini-Pak Size, 14 Fr, 16"	BZK Wipe
42116	16	16"	MTG EZ-Advancer® No-Touch Closed System Catheter Kit, 1500ml Collection Bag, Mini-Pak Size, 16 Fr, 16"	BZK Wipe

MTG EZ-Advancer® Female No-Touch Closed System Catheters

Item #	FR Size	Length	Description	Antiseptic Wipe
Single Closed System Catheters - Pocket-Size				
40412	12	6"	MTG EZ-Advancer® Female No-Touch Closed System Catheter, 700ml Collection Bag, Pocket-Size, 12 Fr, 6"	NA
40414	14	6"	MTG EZ-Advancer® Female No-Touch Closed System Catheter, 700ml Collection Bag, Pocket-Size, 14 Fr, 6"	NA
Closed System Catheters Kits - Pocket-Size				
42412	12	6"	MTG EZ-Advancer® Female No-Touch Closed System Catheter Kit, 700ml Collection Bag, Pocket-Size, 12 Fr, 6"	(2) BZK Wipes
42414	14	6"	MTG EZ-Advancer® Female No-Touch Closed System Catheter Kit, 700ml Collection Bag, Pocket-Size, 14 Fr, 6"	(2) BZK Wipes

MTG EZ-Advancer® Pediatric No-Touch Closed System Catheters

Item #	FR Size	Length	Description	Antiseptic Wipe
Single Closed System Catheters				
30108	8	10"	MTG EZ-Advancer® Pediatric No-Touch Closed System Catheter, 700ml Collection Bag, Pocket-Size, 8 Fr, 10"	NA
30110	10	10"	MTG EZ-Advancer® Pediatric No-Touch Closed System Catheter, 700ml Collection Bag, Pocket-Size, 10 Fr, 10"	NA
Closed System Catheters Kits - Standard Size				
32508*	8	16"	MTG EZ-Advancer® Pediatric No-Touch Closed System Catheter Kit, 1500ml Collection Bag, 8 Fr, 16"	BZK Swabs
32510*	10	16"	MTG EZ-Advancer® Pediatric No-Touch Closed System Catheter Kit, 1500ml Collection Bag, 10 Fr, 16"	BZK Swabs
Closed System Catheters Kits - Mini-Pak Size				
32108	8	10"	MTG EZ-Advancer® Pediatric No-Touch Closed System Catheter Kit, 700ml Collection Bag, Mini-Pak Size, 8 Fr, 10"	BZK Wipe
32110	10	10"	MTG EZ-Advancer® Pediatric No-Touch Closed System Catheter Kit, 700ml Collection Bag, Mini-Pak Size, 10 Fr, 10"	BZK Wipe

*Indicates items that do not include privacy bags.

HCPCS: A4353, Packaging: 100/Case